



# AR

Allround Rollator



For more joy of life:  
Premium quality and fully equipped

## The Saljol Allround Rollator

### Premium quality with many extras

The premium product impresses with its high quality and offers numerous features for your safety, such as over 30 reflectors, non-slip, skin-friendly handles, a practical stick holder, and top-quality brakes. With additional accessories, you can further optimise your rollator.

### Safety you can feel – „hand-hugging“ Cork-TPE handles with armrest

The push handles are made from skin-friendly cork-TPE, which combines the sustainability, lightweight nature, and tactile feel of natural cork with the durability and resistance of plastics. They offer a secure grip thanks to their ergonomic wing shape, preventing your hands from slipping. The shock-absorbing and heat-retaining material reduces pain and promotes the mobility of your hands and fingers, especially for sensitive skin. The non-slip, quick-drying properties ensure secure comfort even in wet conditions. The Cork-TPE armrest makes resting your arms while seated very comfortable. As cork is a natural component of the handle, a typical dark patina may form due to the natural material's reaction with your palms. However, the handles are still easy to clean and remain hygienic.





### **Cork-coloured design elements guide you through the functions of Saljol Rollators**

The design makes intuitive operation easier thanks to cork-coloured, brown components. All adjustable elements are immediately recognisable, allowing quick and secure adjustment to your needs. The well-thought-out colour scheme ensures clarity and maximum comfort in everyday use.



### **Visibility ensures safety**

More than 30 reflectors and reflective strips on the rollator ensure that you remain highly visible in low-light conditions. Whether in traffic or during leisure activities, the reflective elements make you and your rollator more easily noticeable to other road users.



### **Optional one-hand braking system – brake with one or both hands for maximum control at no extra charge (on request)**

The conversion to a one-hand simultaneous brake allows you to fully brake the rollator with just one hand. This braking system is particularly helpful for those experiencing age-related loss of strength or unilateral limitations, such as after a stroke, with Parkinson's, or fractures, and enables you to operate both brakes with minimal effort. The necessary longer cable for the brake system will be sent to you free of charge upon request.



### **Compact and foldable – ideal for easy transport**

The rollator can be folded to a width of just 21 cm with a single movement. When folded, it can be securely locked with a hook, making it easy to load and store compactly and stably. Thanks to its lightweight design and narrow dimensions, it can be easily transported in a car or while travelling.



### **Highest quality in every detail – you'll feel it!**

The smooth-running hand brakes allow for effortless braking and are easy to operate, even with smaller hands. The ergonomic shape and cork-TPE covering provide a secure and comfortable grip. When you engage the parking brake, a clear clicking sound lets you know that the rollator is securely in place – ideal for sitting down or filling the bag.



#### **Extra large, shock-absorbing and grippy – wheels for every situation**

The tyres with a particularly grippy tread provide stability and comfort on uneven surfaces such as cobblestones, field, or forest paths. In addition, the 250 mm large all-terrain front wheels prevent the wheels from getting blocked by small obstacles such as roots or potholes. This reduces the risk of falling and allows for safe movement – ideal if your sense of balance is impaired or your attention is reduced. The shock-absorbing PU solid rubber material functions almost like an air tyre, but without the risk of punctures. It also protects your hands, elbows, and shoulders, providing noticeable relief for those with rheumatic or neurological conditions.



#### **Overcome obstacles safely and effortlessly with non-slip tipping aids**

The non-slip, rubberised tipping aids are positioned at the rear of the frame, near the back wheels. They make it easier to navigate curbs, thresholds, and uneven surfaces. The flat step surface provides a secure grip, reduces the risk of falls, and offers added stability.



**Secure hold for your walking stick – always within reach**

The universal cane holder adapts flexibly to different mobility aids and walking sticks. With a unique locking system, the stick is securely fixed and always within reach. A clicking sound lets you know that the stick is safely attached.



**Your support wherever you go – the lightweight and foldable walking stick**

The walking stick is incredibly lightweight at just 278 g, yet remarkably durable. Made from high-quality carbon fibre, it impresses with its top quality and elegant look. It is also height-adjustable, foldable and always stable – the perfect combination of functionality and elegance. It is available directly at [www.saljol.de/en](http://www.saljol.de/en), through our service team, or at a medical supply store near you.



**The right product for you – available for every need**

Saljol Premium Rollators are available in two sizes, elegant colours, and additionally with a wider seat (52 cm) and higher weight capacity (180 kg). This gives you the option to choose the solution that best suits your needs and personal taste. Whether for more comfort, functionality, or style – you will find the perfect combination to meet your requirements.



**Retrofittable grinding brake – your rollator can be tailored to your needs**

Safety and control while walking are your top priorities? The retrofittable grinding brake (accessory for your rollator) provides exactly that. The rolling resistance of the rear wheels can be individually adjusted, gently braking without locking (ABS), and preventing your rollator from rolling faster than you want – especially helpful on slopes. The brake supports you optimally in cases of walking instability, can counteract incorrect directional drift, and thus reduces the risk of falling. Once properly adjusted, you can activate or deactivate it with your foot.





German Engineering from Bavaria.

**German Design**



**Great luxury that remains affordable**

Saljol Rollators offer the perfect mix of functionality, quality, and a touch of luxury. German engineering combines high-quality features with comfort to make your everyday life easier – all at a fair overall price.



**Comprehensive support – for your mobility and safety**

Our physiotherapist, Katharina, supports you in all areas where expert physiotherapeutic knowledge is required. Our helpful customer service is always at your side – whether in the medical supply store or directly with us. If you have any questions or uncertainties regarding your Saljol Allround Rollator, we provide the support you need to get the most out of your rollator.



**Training with your rollator**

To avoid medical side effects when using our rollators, we invented Rollastic. It's a tailored exercise programme for rollator users that, when performed regularly, can help maintain mobility and strength. Special videos are provided with the rollator. For your safety in road traffic, we developed the safety training „Rollator Mobility“. Helpful tips and tricks can be found at [www.rollatormobility.com](http://www.rollatormobility.com).



**Peter: „I've got my mobility back!“**

„Thanks to the Allround Rollator, I regained much of my mobility after my stroke. The lack of mobility in my left hand is a problem for me, and I struggle immensely with fine motor skills. The one-hand brake on the Allround Rollator helps me enormously. After the stroke, I came home completely helpless and thought that travelling would no longer be possible. Wrong! There's always a way! You just have to give it a try.“



| Allround Rollator                | AR54       | AR62       | AR62W      |
|----------------------------------|------------|------------|------------|
| Recommended body height          | 145–170 cm | 160–200 cm | 160–200 cm |
| Maximum user weight              | 150 kg     | 150 kg     | 180 kg     |
| Total weight without accessories | 6.3 kg     | 6.4 kg     | 6.4 kg     |
| Maximum load of net bag          | 5 kg       | 5 kg       | 5 kg       |
| Width open/closed                | 61/21 cm   | 61/21 cm   | 68/21 cm   |
| Seat width                       | 46 cm      | 46 cm      | 52 cm      |
| Seat height                      | 54 cm      | 62 cm      | 62 cm      |
| Handle height (adjustable)       | 76–88 cm   | 81–103 cm  | 81–100 cm  |


### Allround Rollator accessories



Back belt



LED light

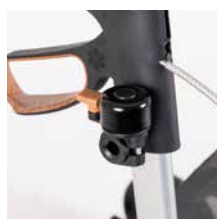

Bell



Stick holder



Grinding brake



Bamboo tray



Firm seat



Travel bag



Umbrella  
caramel



Umbrella  
black



Walking stick




Brake cable for  
the conversion to  
one-hand brake

### Frame colour



Pearl Grey

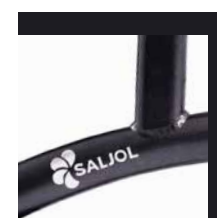

Magic Black



Forest Green



Barolo Red



Saljol GmbH | [www.saljol.de/en](http://www.saljol.de/en)  
Frauenstraße 32 | D-82216 Maisach

